


Flavor by Japanese Standards of Food Additives

Data No. GA230-0591

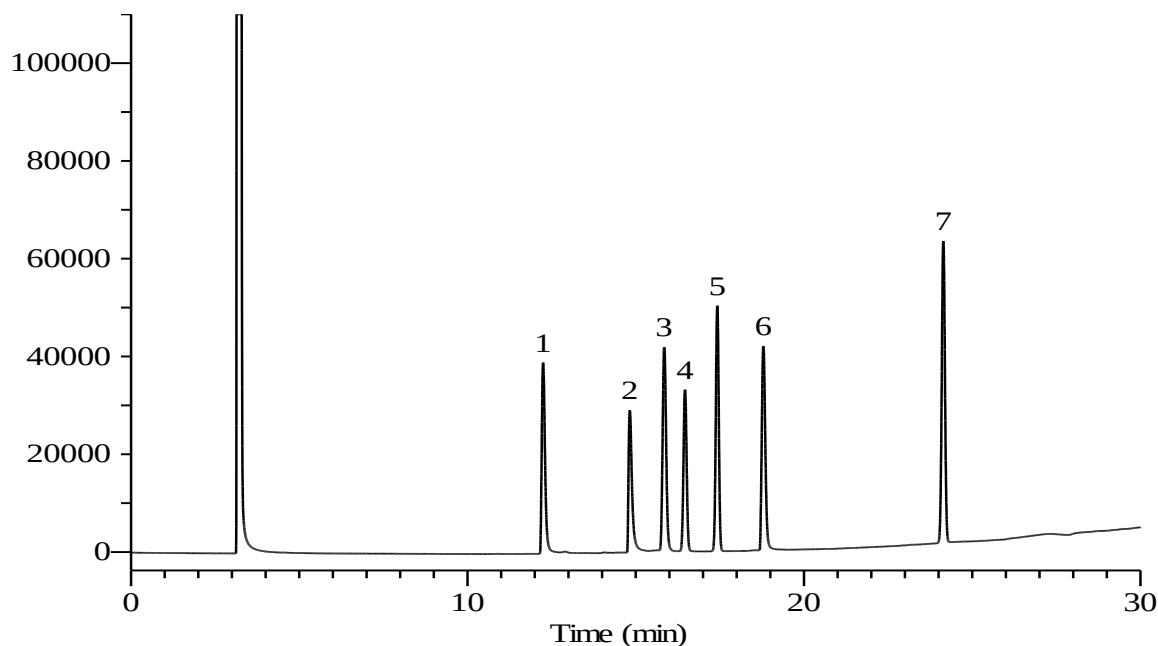

Conditions

System : GC-FID
Column : InertCap 1
0.32 mm I.D. x 30 m df = 5.00 µm
Col. Cat. No. : 1010-11249
Col. Temp. : 50 °C(5 min hold) - 5 °C/min - 230 °C
Carrier Gas : He 50 kPa
Injection : Split flow 90 mL/min
175 °C
Detection : FID Range 10⁴
250 °C
Sample Size : 1000 µg/mL in Methanol
1 µL

Analyte

1. Butylamine
2. Pyrrolidine
3. Isopentylamine
4. Pyrazine
5. Pyrrole
6. Piperidine
7. 2,6-Dimethylpyridine