

InertSearch™ for LC


Inertsil® Applications

Analysis of Vitamin C (Inertsil HILIC)

Data No. LA632-0696

Analyte:

1. Erythorbic acid (5 mg/L)
2. Ascorbic acid (5 mg/L)

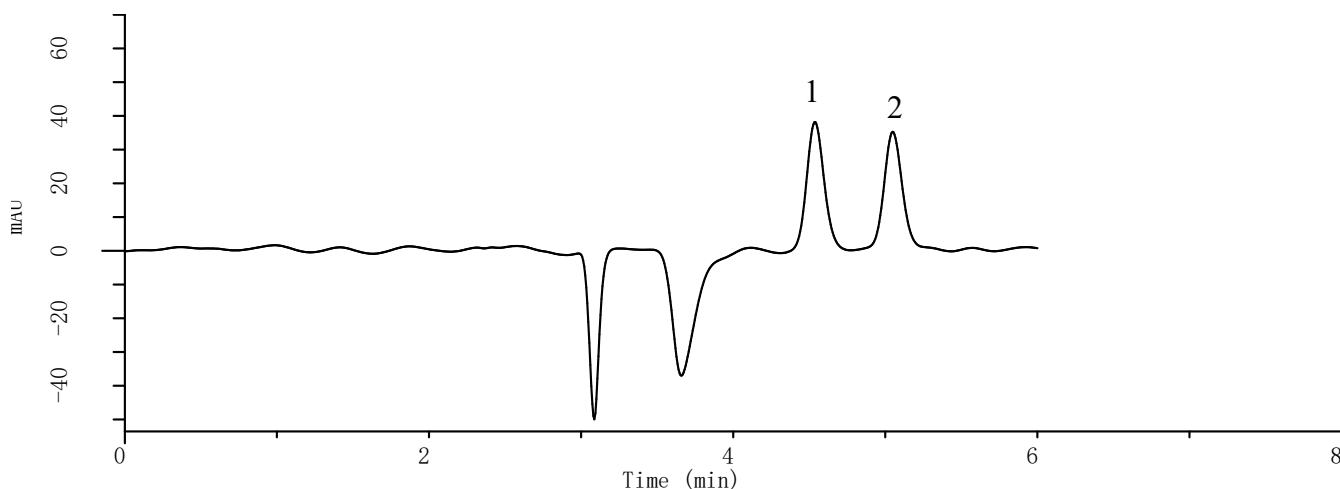


Conditions


System : GL-7400 HPLC system
Column : Inertsil HILIC (5 μ m, 250 x 4.6 mm I.D.)
Column Cat. No. : 5020-07736
Eluent : A) CH₃CN
 B) H₂O
 C) CH₃COOH
 A/B/C = 87/11/2, v/v/v
Flow rate : 1.0 mL/min
Col. Temp. : 40 °C
Detection : UV 243 nm (GL-7452A PDA Detector)
Injection Vol. : 20 μ L
Sample : Vitamin C

Analysis of Vitamin C (Inertsil HILIC)

Data No. LA632-0696

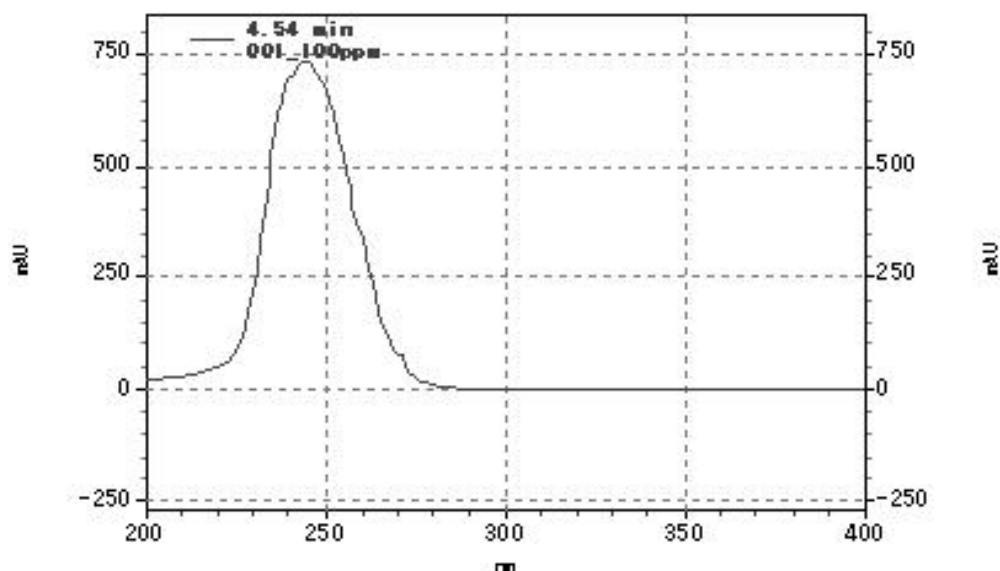

UV spectrum of Erythorbic acid (ErA)

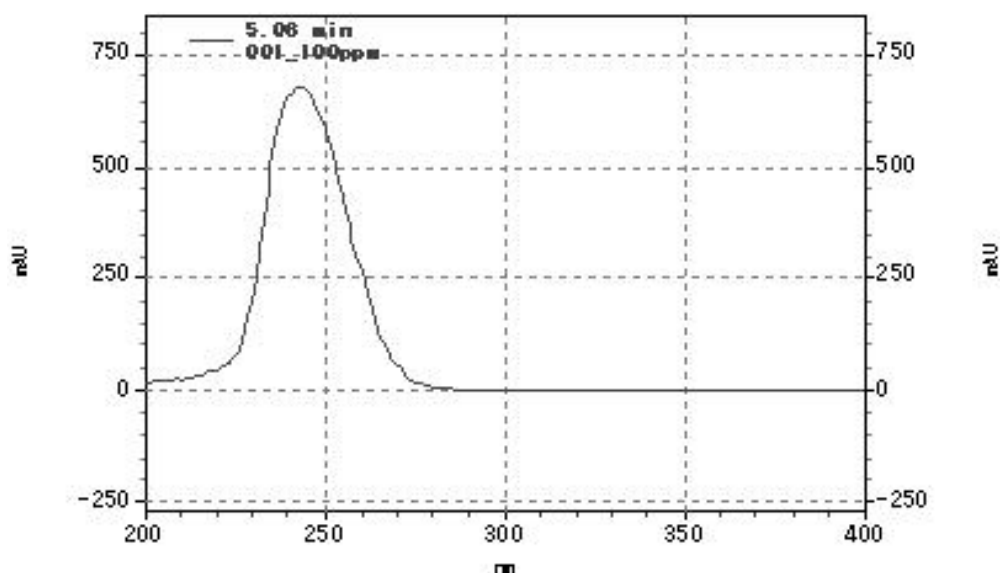

UV spectrum of Ascorbic acid (AsA)