


Analysis of Nucleotide

Data No. LA683-0848

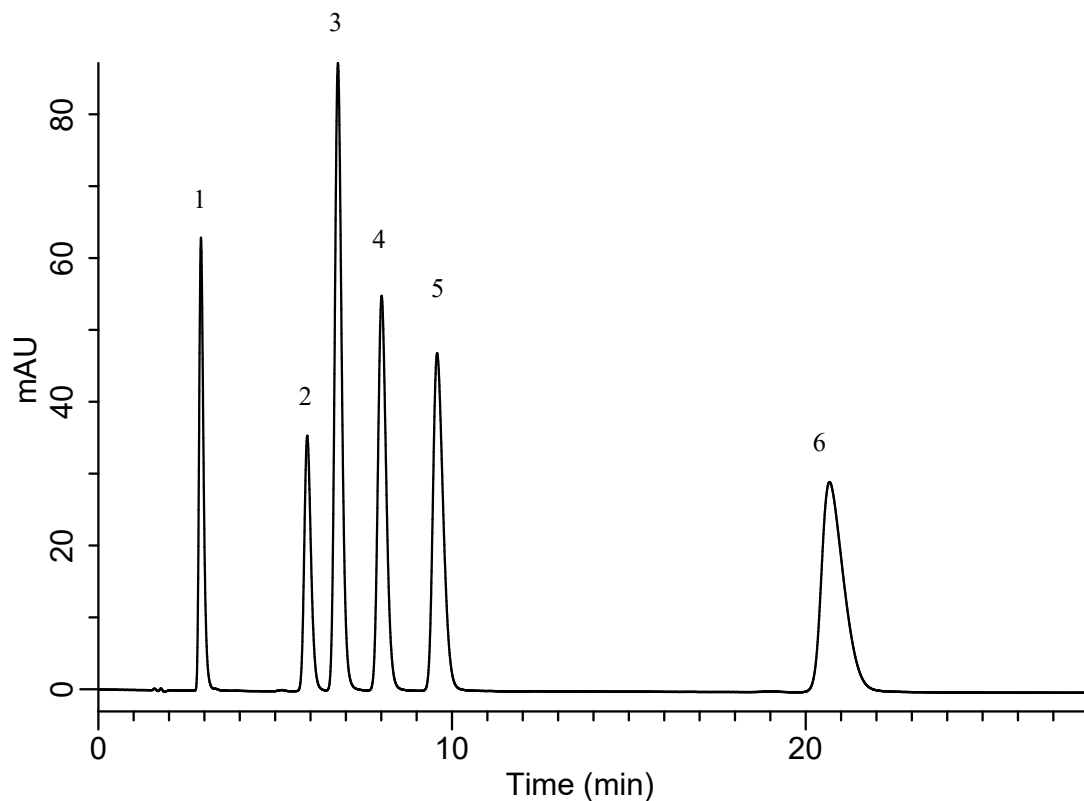

Conditions

System : GL-7400 HPLC System
Column : Inertsil ODS-3 (5 μ m, 150 x 3.0 mm I.D.)
Column Cat. No. : 5020-01751
Eluent : 5 mM KH_2PO_4 (pH 2.5, H_3PO_4)
Flow Rate : 0.4 mL/min
Col. Temp. : 40 $^{\circ}\text{C}$
Detection : UV 260 nm (GL-7450 UV Detector)
Injection Vol. : 10 μ L

Analyte

Nucleotides
1. 5'-CMP (38 $\mu\text{g/mL}$)
2. 5'-UMP (38 $\mu\text{g/mL}$)
3. 5'-AMP (38 $\mu\text{g/mL}$)
4. 5'-GMP (38 $\mu\text{g/mL}$)
5. 5'-IMP (38 $\mu\text{g/mL}$)
6. 5'-TMP (38 $\mu\text{g/mL}$)